


TABLET / CAPSULE PRINTING MACHINE

“Print to BRAND”



Engineering Solutions For Healthcare

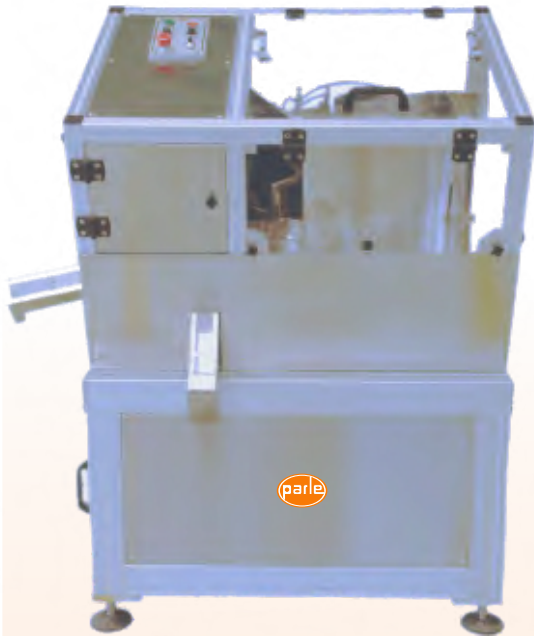

PPT - 251




PPT - 321

PARLE GLOBAL TECHNOLOGIES PVT. LTD.

www.parleglobal.com | info@parleglobal.com

Salient Features


- » Automatic printing Machine for Tablet Capsule and soft gelatin
- » Minimum change parts, equipped with feed discs and gravure rollers
- » Can do linear printing on any shape (round, triangular, oval, spherical or hexagonal)
- » Special designed discharge chute arrangement for sorting the tablets printed during the setting.
- » Interlocks for low air pressure, machine stop.
- » Total sound free machine with silent drive transmission.
- » Integral variable speed controls.
- » Operator friendly.
- » All contact parts in approved food grade MOC, confirming cGMP standards with safety standards.
- » Fine adjustable type printing head, for aligning the print on tablet center.
- » Fine adjustable pressure roller setting, giving constant pressure transfer on tablet/capsule to be printed.

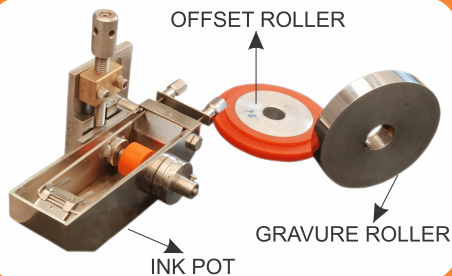

AIR FILTER & REGULATOR



MACHINE DRIVE

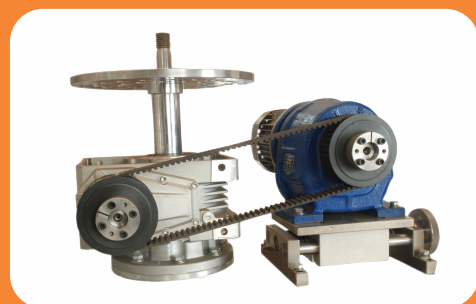

MAIN DRIVE ASSEMBLY

Application



FOR PRINTING ON
TABLETS & CAPSULES

System Basic Requirement

- All tablets surfaces need to be either film / enteric / sugar coated, polished and dust free.
- Soft gel capsules should be oil free.
- Food grade ink and approved IPA required for printing process.

Technical Specification

Description / Model		PPT-25 I	PPT-32 I
Output	Hard / Soft	55,000 per hour (#0)	75,000 per hour (#0)
	Gelatin Capsules		
	Caplets / Tablets	65,000 per hour Ø 11mm	85,000 per hour Ø 11 mm
Power Supply		230 V, 1ph, 50 / 60 Hz	230 V, 1ph, 50 / 60 Hz
Power Consumption		0.4 KW, 0.5 H.P.	0.4 KW, 0.5 H.P.
Air Consumption Pressure		40 lit. / min. / 6 kg/cm ²	40 lit. / min. / 6 kg/cm ²
Net weight		480 Kg.	650 Kg
Machine Dimensions		1025 (L) x 1100 (B) x 1450(H) mm	1210(L) x 1270 (B) x 1300 (H) mm
Gross weight		525 Kg.	725 Kg.
Gross Dimensions		1525 (L) x 1350 (B) x 1780 (H) mm	1625 (L) x 1725 (B) x 1825 (H) mm

Our Representatives in USA:

Daiichi Jitsugyo (America), Inc.

939 AEC Drive, Wood Dale, IL 60191

Phone : 630-776-2459

Email: Rnelson@djausa.com

Web: www.dja-pharma.com • www.parleglobalusa.com

